



Streamline workflow in your office.

The Welch Allyn CP 200™ Electrocardiograph

WelchAllyn®

Advancing Frontline Care™

Welch Allyn CP 200 Electrocardiograph— An easy-to-use, full featured device with

Welch Allyn CP 200™

- > A non-interpretive 12-lead resting ECG that provides accurate measurements
- > Optional MEANS ECG interpretation software offers a second opinion
- > Full alphanumeric keyboard and LCD screen allow you to enter patient data quickly
- > High-resolution thermal printer generates standard full-sized reports in the format you need
- > Help feature lets you access patient education tools as well as answer system operation questions
- > Uses standard SD cards for virtually unlimited ECG data storage or to transfer data to the Welch Allyn CardioPerfect™ Workstation Software for analysis and/or transfer to an EMR
- > Colour ECG waveform preview screen makes it easier to read critical information
- > Dedicated function keys speed access to common/typical actions
- > Store up to 50 ECGs and 50 spirometry tests for later review and editing of patient data and interpretation, or add SD cards for virtually unlimited data storage
- > Transfer patient data to the Welch Allyn CardioPerfect™ Workstation Software and your EMR via SD cards or via direct connection
- > Use Schedule Test mode to programme all patient information in at the beginning of the day so you can use your time more efficiently
- > Administrative functions with password log-on controls to support your program for HIPAA compliance
- > Optional spirometry testing and optional MEANS interpretation ECG software



Colour LCD Screen
Large colour LCD screen allows you to preview information on screen as well as review tests before printing



Dedicated Function Keys
Save time by using "hot" keys to perform the most common functions without looking through menus

Components/Accessories of the CP 200

Optional ECG Cart

- > Shelf directly under the unit for paper or manuals
- > Attractive, ergonomic design
- > Hospital-grade wheels also available

Cable Arm and Shelf Option



Digital ECG Patient Cable

- > Lead wires match their anatomical positions, speeding patient hookup while reducing tangling
- > Individually replaceable lead wires are more cost-effective



Data Storage on SD Cards

- > Virtually unlimited storage of ECG and spirometry tests
- > Transfers data to the Welch Allyn CardioPerfect™ Workstation Software for easy retrieval, review, and storage



advanced options to help streamline workflow.



Direct Connect Connectivity

Connect to the Welch Allyn CardioPerfect™ Workstation Software and your EMR easily through a cable interface



Welch Allyn CP 200™ Combo – ECG with optional spirometry capability

The Welch Allyn CP 200 with optional spirometer utilizes one platform for two different testing applications. You'll enjoy all the benefits of a full-featured ECG and spirometer with data organised in the same fashion to make it easier to manage.

- > Performs FVC, SVC and pre- and post-BD testing with automatic interpretation and reversibility analysis
- > Improve test quality and patient coaching with our on-screen quality messages and variability checks
- > Educate your patients with our printed help tool summarising their spirometry results
- > Customise report formats and screen displays to fit your needs and workflow
- > Spirometry tests meet all industry standards including ATS, OSHA, NIOSH, NLHEP and Social Security
- > Perform single or multiple flow calibration with automatic printing to save time
- > Automatic lung age calculation to help with smoking cessation programmes



Use the Welch Allyn CardioPerfect™ Workstation Software to manage your ECG data

Easy-to-use, intuitive software that is EMR ready

The Welch Allyn CP 100™ and CP 200™ transfer information to the Welch Allyn CardioPerfect™ Workstation Software, which provides paperless efficiency as well as connectivity to many leading EMR systems. This software gives you the ability to store, access, review, compare and edit diagnostic test data from any computer on your network.


- > Retrieve records with the touch of a key, and review and edit right on the screen
- > Store thousands of records with ease
- > Access information from your exam rooms, your office, or at home
- > Fax or e-mail reports directly for consultant input
- > Integrate with other patient information systems



Welch Allyn CP 100™
Transfer information via SD card



Welch Allyn CP 200™
Transfer information via direct connection or SD card

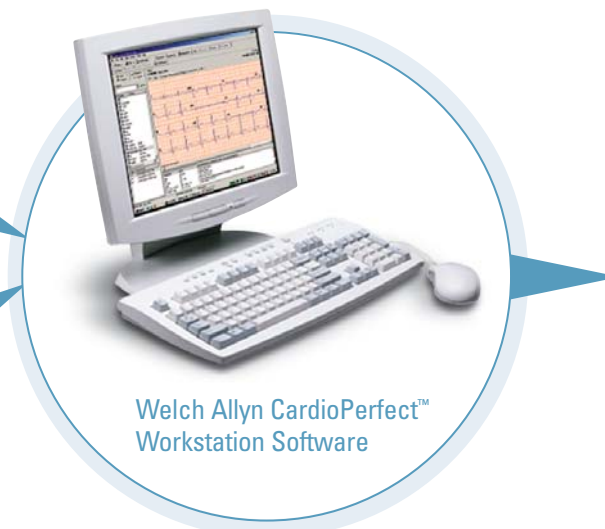

Welch Allyn CardioPerfect™ Workstation Software



Transfer test data and patient information to individual patient records in your EMR for analysis and permanent storage

Ordering Information

United Kingdom	CP 200	CP2-2E4	CP 200 with IEC Patient Cable, English, UK Power Cord
	CP 200i	CP2A-2E4	CP 200 with IEC Patient Cable, English, UK Power Cord, Interpretation
Europe/Middle East	CP 200	CP2-2E2	CP 200 with IEC Patient Cable, English, Eur Power Cord
	CP 200i	CP2A-2E2	CP 200 with IEC Patient Cable, English, Eur Power Cord, Interpretation
South Africa	CP 200	CP2-2E9	CP 200 with IEC Patient Cable, English, South African Power Cord
	CP 200i	CP2A-2E9	CP 200 with IEC Patient Cable, English, South African Power Cord, Interpretation

ECG Accessories

45008-0000	Resting Tab Electrodes (Box of 1000)
58581-0000	Resting Tab Electrode Adaptors (Set of 10)
94018-0000	Paper, 200 sheets/pad, 5 pads/case
401393	Office Cart
401394	Hospital Cart
401161	Cable Arm and Shelf Option
400294	ECG Patient cable, IEC
401122	Leadwire Set (10), IEC, Banana
401124	Leadwire Set (10), IEC, Pinch
100660	ECG Replacement Battery Assembly
701586	Dust Cover
100623	Interpretation Upgrade
71300-0000	ECG Hookup Wall Chart

Spirometry Accessories

100400	Spirometry Upgrade Kit, with 3 L Calibration Syringe
100680	Nose Clip (25 Pieces)
703480	Calibration Syringe (3 L)
703418	Disposable Flow Transducers, 25 each
703419	Disposable Flow Transducers, 100 each
703415	Pressure Tubing, 2 meters
703552	Replacement Sensor, RS232
71038-3000	Spirometry Reference Chart – English
703337	Effort Acceptability Poster – English

Welch Allyn U.K.

Cublington Road
Aston Abbotts
Buckinghamshire HP22 4ND
England
Telephone: +44 (0) 207 365 6780
Fax: +44 (0) 207 365 9694

Welch Allyn South Africa (Pty) Ltd

394 Surrey Ave, Ferndale, Randburg 2194
South Africa
PO Box 809, Fontainebleau, South Africa, 2032
Telephone: +27 (0) 11 777 7555
Fax: +27 (0) 11 777 7556

Welch Allyn Europe

Navan Business Park
Dublin Road
Navan, County Meath
Ireland
Telephone: +353 46 906 7790
Fax: +353 46 906 7754